

Chapter 20

In Vitro Dissection of Protein Translocation into the Mammalian Endoplasmic Reticulum

Ajay Sharma, Malaiyalam Mariappan, Suhila Appathurai,
and Ramanujan S. Hegde

Abstract

In eukaryotic cells, roughly one-fourth of all mRNAs code for secretory and membrane proteins. This class of proteins must first be segregated to the endoplasmic reticulum, where they are either translocated into the lumen or inserted into the lipid bilayer. The study of these processes has long relied on their successful reconstitution in cell-free systems. The high manipulability of such in vitro systems has allowed the identification of key machinery, elucidation of their functional roles in translocation, and dissection of their mechanisms of action. Here, we provide the basic methodology for (i) setting up robust mammalian-based in vitro translation and translocation systems, (ii) assays for protein translocation, insertion, and topology, and (iii) methods to solubilize, fractionate, and reconstitute ER membranes. Variations of these methods should be applicable not only to forward protein translocation systems but also for dissecting other poorly understood membrane-associated processes such as retrotranslocation.

Key words: in vitro translation, microsomes, membrane proteins, reconstitution, proteoliposomes, protease protection, protein topology.

1. Introduction

The endoplasmic reticulum (ER) is the major site for the biosynthesis, maturation, quality control, and degradation of secretory and membrane proteins. Each of these basic processes employs multiple distinct pathways that operate in parallel to provide the cell considerable flexibility in handling the tremendous diversity of proteins that transit through the ER. A major goal of cell biology has long been to identify and dissect the mechanism of action of the machinery that define these pathways of secretory and membrane protein metabolism.

One of the principal approaches that has been applied to this problem is the reconstitution of key pathways or sub-reactions in a cell-free system. The tremendous manipulability of such systems affords a direct window into the biochemical and mechanistic dissection of any of these pathways. By understanding the basic features of the standard cell-free translation and translocation system, it can be sensibly customized in various ways to address any of a number of processes that occur at or within the ER and cytosol. A few of the major ER-associated processes that are amenable to dissection using this or similar *in vitro* systems are shown in **Fig. 20.1**.

In co-translational translocation (**Fig. 20.1A**), secretory and membrane proteins are recognized as they are being synthesized. The signal recognition particle (SRP) binds to a hydrophobic

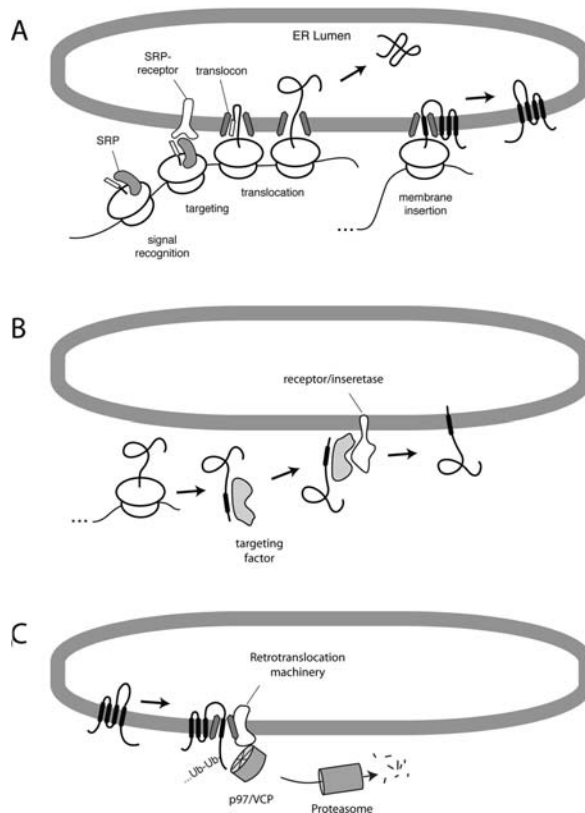


Fig. 20.1. Examples of ER-associated pathways amenable to *in vitro* reconstitution. **(A)** The SRP-dependent co-translational translocation pathway. **(B)** A post-translational translocation pathway for tail-anchored membrane protein insertion. **(C)** Post-translational pathway of membrane protein metabolism involving ubiquitination, retrotranslocation, and proteasomal degradation. In each of these instances, the substrate is synthesized *in vitro*, making it the only protein that becomes radiolabeled. The other components of the system can be manipulated to analyze the requirements for substrate translocation, maturation, degradation, etc.

domain (either a signal sequence or transmembrane domain) in the nascent polypeptide as it emerges from the ribosome. This complex of ribosome–nascent chain–SRP is then targeted to the membrane (via the SRP receptor) and transferred to a translocon whose central component is the Sec61 complex. Secretory proteins are translocated through the Sec61 complex, while membrane proteins are laterally released by the Sec61 complex into the membrane bilayer. Because everything occurs co-translationally, the study of these events typically depends on translating the protein of interest in a cytosolic extract *in the presence of* a source of ER membrane [typically rough microsomes (RM) isolated from canine pancreas (1) or other tissues/cells (2, 3)]. If the ER membranes are added to the reaction after translation is completed, translocation will not occur.

In post-translational translocation pathways (**Fig. 20.1B**), the protein remains competent for translocation even after it has been fully synthesized and released from the ribosome. Such pathways can be employed by a subset of secretory proteins, some small proteins, and certain types of membrane proteins. These pathways, particularly in higher eukaryotes, are not nearly as well understood as the SRP-dependent co-translational translocation pathway. For example, a novel and well-conserved pathway for tail-anchored membrane proteins was only recently discovered (4, 5) and whose full complement of machinery remains to be clarified. From a practical standpoint, these post-translational pathways are in many ways easier to study because the translation reaction can be uncoupled from the translocation reaction. Thus, proteins can be translated, after which various manipulations can be applied (e.g., removal of energy, addition of inhibitors, change in conditions) before initiating the translocation reaction by adding a source of ER membranes. This provides greater flexibility than co-translational reactions where conditions must be maintained within the narrow range that is compatible with efficient protein synthesis.

And finally, reactions after both synthesis and translocation have been completed can also be studied using these same in vitro systems (**Fig. 20.1C**). Examples of such processes include maturation events in the ER lumen (6), quality control of misfolded proteins, retrotranslocation, ubiquitination, and degradation (7). Again, many of these pathways are still relatively poorly understood, particularly from a mechanistic point of view. As with post-translational translocation, these events can often be uncoupled from protein synthesis (and in some cases, even translocation), allowing experimental flexibility.

The study of all of these processes in vitro depends on three basic tools. An in vitro translation (IVT) system, clear and definitive assays for translocation and topology, and, in the case of events occurring at the ER, methods to manipulate the

composition of the membrane. These basic tools can be applied in a wide range of ways. The IVT system allows one to produce in a physiologic system a protein radiolabeled with very high specific activity that can be followed. By scaling up these reactions, biochemical amounts can also be generated to identify interacting partners (4). Highly specific assays for translocation can be used to identify protein or lipid requirements (4, 8, 9) and analyze the action of small molecule inhibitors (10). Manipulation of the membrane provides access to the requirements at this usually inaccessible compartment (8, 9, 11). And the ability to isolate the membrane after insertion provides the ability to study subsequent events (such as degradation) in isolation (7).

2. Materials

2.1. Preparation of the Transcription Mix (T1)

For general advice regarding these reagents, see the following:

1. 1 M HEPES, pH 7.6: Prepare a solution of 1 M HEPES (free acid), titrated with 0.45 M NaOH. This will be the correct pH when diluted in the buffers below. Filter and store at 4°C (*see* **Notes 1, 2, 8, and 9**).
2. 2 M MgCl₂, store at RT or 4°C.
3. 100 mM spermidine (Sigma); very hygroscopic. Dissolve 145 mg spermidine in 10 mL water. Freeze in nitrogen and store at -20°C.
4. 1 M DTT (1,4-dithiothreitol; Roche). Dissolve 1.54 g DTT in 10 mL water. Aliquot and freeze in nitrogen. Store at -80°C. Do not freeze-thaw more than twice to prevent oxidation.
5. 10X NTPs. 5 mM each of ATP, UTP, and CTP and 1 mM of GTP in water. Adjust to pH ~7 with NaOH as needed. Aliquot and freeze in nitrogen. Store at -80°C. Do not freeze-thaw more than five times.
6. 10X Cap: 7-methyl diguanosine triphosphate cap structure analog (New England Biolabs). Each vial contains 25 A₂₆₀ units. Add 300 µL water directly to the vial (to make ~5 mM solution), mix well to dissolve, aliquot, and freeze in nitrogen. Store at -80°C. Do not freeze-thaw more than five times. (*see* **Note 3**).

2.2. Preparation of the Translation Mix (T2)

1. Crude rabbit reticulocyte lysate (RRL). This can be prepared in-house (12), although very few labs currently do this due to practical limitations. Purchase from Green Hectares. They offer two products. We buy the more dilute material which they say is for 'purification of factors.' The

other product 'for in vitro translation' is somewhat more concentrated (and expensive), but we have not found this to be a significant advantage. The more dilute product typically comes in rather large aliquots (~50 mL), which should be stored at -80°C . Do not freeze-thaw more than twice. The first time you thaw a 50 mL aliquot, make 10 mL aliquots in 15 mL polypropylene tubes. To one of these 10 mL aliquots, add hemin, treat with micrococcal nuclease, and further sub-aliquot it (*see Section 3.2.*). The other 10 mL aliquots can be frozen directly in nitrogen and stored at -80 until needed. Crude lysate is stable for up to 10 years. Nuclease-treated RRL may be less stable, which is why we do not nuclease everything at once.

2. Micrococcal nuclease (Calbiochem; 15,000 units/vial). Dissolve in 1 mL 50 mM HEPES buffer pH 7.4 (i.e., 15 units per μL); aliquot and freeze in liquid nitrogen and store at -20°C .
3. Hemin (Sigma). Prepare 10 mL of a 100 μM stock by mixing the following *in this exact order* (to avoid problems with precipitation): 6.5 mg hemin, 250 μL 1 N KOH, 500 μL 200 mM Tris, pH 8.0, 8.9 mL ethylene glycol, 190 μL 1 N HCl, 50 μL water. Mix by vortexing and store at -20°C ; will not freeze due to ethylene glycol.
4. 100 mM CaCl_2 : Dissolve in water. Keep at 4°C .
5. 200 mM EGTA. Prepare as follows: To 760 mg EGTA powder, add ~7 or 8 mL water and vortex to resuspend the powder (it will not go into solution). Add 950 μL 5 N NaOH (at which point the EGTA will go into solution). Adjust to 10 mL with additional water. Store at 4°C .
6. Calf liver tRNA (Novagen). This is typically supplied as a 10 mg/mL stock. Freeze in nitrogen and store at -80°C . Stable to multiple freeze-thaws as long as it is kept on ice when thawed and frozen in liquid nitrogen immediately after use (*see Note 4*).
7. 1.2 M creatine phosphate (Roche). Dissolve in water, freeze in nitrogen, and store at -80°C .
8. 20 mg/mL creatine kinase (Roche). Dissolve 100 mg in 5 mL of 10 mM HEPES, pH 7.5, 50% glycerol. Store at -20°C .
9. 5 M KOAc stock (for each 500 mL add 2 mL of 12 N HCl to bring pH to ~7). Store at RT.
10. Amino acid stocks and mixes – each one is made up individually as a 20 mM stock, and these are mixed to prepare an amino acid mix of 19 amino acids (1 mM each) without methionine (which will be supplied as a ^{35}S -methionine to

label translated proteins). Prepare stock solutions of each amino acid (purchased as powders from Sigma) at 20 mM in either 0.01 N HCl (Trp, Val, Ile, Asn, Phe, Asp, Glu, and Lys), 0.1 N HCl (Tyr), or water (the remaining ones). To make the 19 amino acid mix (1 mM each) simply mix 0.5 mL of each of the 19 with 0.5 mL of H₂O. Mixes lacking different amino acids for labeling with residues other than methionine can also be made. Freeze in aliquots and store both the stock amino acids and the mixes at -80°C (*see* **Note 5**).

11. 10X Emix: for 1 mL, mix 305 μL water, 400 μL 19 amino acid mix (1 mM each; step 10 above), 10 μL 1 M HEPES (not pH adjusted), 1.87 μL 8 N KOH, 83.3 μL 1.2 M creatine phosphate (step 7 above), 100 μL 0.1 M ATP, 100 μL 0.1 M GTP. Freeze in nitrogen and store at -80°C .
12. CB20X: 627.5 μL water, 240 μL 1 M HEPES (not pH adjusted), 22.5 μL 8 N KOH, 100 μL 5 M KOAc pH 7 (step 9 above), 10 μL 2 M MgCl₂. Freeze in nitrogen and store at -80°C (*see* **Note 10**).

2.3. Linked Transcription and Translation

1. Reagents for PCR amplification of the desired cDNA.
2. Qiagen PCR purification kit.
3. RNAsin (Promega). *Do not freeze in nitrogen*; keep at -20°C . This comes in 50% glycerol. Stable for at least 2 years.
4. T7 or SP6 polymerase (New England Biolabs). *Do not freeze in nitrogen*; keep at -20°C . This comes in 50% glycerol. Stable for at least 2 years.
5. ³⁵S-methionine (from PerkinElmer; 1,000 Ci/mmol, in aqueous solution). Store in aliquots of 100 μL or less at -80°C . Suitable precautions should be taken when working with radioisotopes.
6. Canine pancreatic rough microsomes (RM). RMs can be obtained in small amounts (50–200 μL) from commercial sources (e.g., Promega). However, for large amounts as needed for fractionation and purification of membrane proteins, purchasing RMs becomes prohibitively expensive (*see* **Notes 6 and 11**).

2.4. Assays for Translocation and Topology

1. Physiologic salt buffer (PSB): 100 mM KOAc, 2 mM Mg(OAc)₂, 50 mM HEPES, pH 7.4. It is often convenient to prepare a 10X PSB stock, which is diluted with various other components (e.g., sucrose) as needed.
2. Proteinase K (PK; Roche): 10 mg/mL dissolved in 20 mM HEPES, pH 7.4. Prepare single-use 20 μL aliquots, freeze in nitrogen, and store at -80 . After thawing and use, discard

remainder. We usually prepare ~100 aliquots at a time (*see Note 7*).

3. 10% Triton X-100 solution. Store at 4°C.
4. PK-kill buffer (PKB): 1% SDS, 0.1 M Tris, pH 8.0.

2.5. Reconstitution of Membrane Proteins into Proteoliposomes

1. Bovine liver phosphatidylcholine (PC, supplied in organic solvent, usually chloroform or chloroform:methanol; from Avanti Polar Lipids) (*see Note 12*).
2. Bovine liver phosphatidylethanolamine (PE, supplied in organic solvent, usually chloroform or chloroform:methanol; from Avanti Polar Lipids).
3. Lissamine–rhodamine–dipalmitoyl–PE (DPPE, supplied in organic solvent, usually chloroform or chloroform:methanol; from Avanti Polar Lipids).
4. Lipid hydration buffer: 50 mM HEPES, pH 7.4, 10 mM DTT (added from a freshly prepared 1 M DTT stock solution).
5. Bio-Beads SM2 (Bio-Rad). Two grades are sold. We buy standard grade product (less expensive) as it seems to behave identically to the more expensive ‘biotechnology’ grade Bio-Beads.
6. RM (*see Note 6*).
7. 10% DeoxyBigChap (DBC; Calbiochem). Add ~9.5 mL water to a 1 g vial of DBC to bring the total volume to 10 mL. Dissolve by gentle agitation at room temperature. Store at 4°C for short-term use (6 months or less); otherwise freeze and store at either –20 or –80°C (*see Note 13*).
8. Pre-extraction buffer: 50 mM HEPES, pH 7.4, 250 mM sucrose, 2 mM MgCl₂, 0.2% DBC.
9. Extraction buffer: 400 mM KOAc, 5 mM MgCl₂, 50 mM HEPES, pH 7.4, 15% glycerol, 1 mM DTT (added just before use from a freshly prepared 1 M stock solution).

3. Methods

In our lab, we have simplified the transcription and translation reactions to essentially be in-house generated ‘kits’ that are easy to use. The kit components consist of a T1 mix (for transcription) and a RRL-based T2 mix (for translation). Because a single (meticulous and responsible) person can be charged with preparing and maintaining aliquots of T1 and T2, the experiments become highly reproducible and sufficiently straightforward that

they are amenable to even the most inexperienced lab members. Furthermore, the use of this approach is far more economical and allows for much larger scale applications of in vitro translation reactions than is possible from commercial kits. This scalability actually makes it reasonable to purify the translated proteins [and identify interesting co-associating factors (4)], something that is typically not considered an option with mammalian in vitro systems.

The stock solutions in Steps 1 and 2 of **Section 2.1** are all quite stable, and we have one primary person responsible for preparing and maintaining them. They are not used for any other purpose and are therefore kept in a separate place for use only in preparing transcription and translation reagents. This is to maintain reliability and prevent RNase contamination. The T1 and T2 mixes are also quite stable, but less so than individual reagents. We therefore make it from the stock solutions from time to time. If you are doing lots of in vitro transcription and translation reactions, you can make 7.6 mL and 17.1 mL T1 and T2 mixes, respectively, at a time (as we usually do). Otherwise, 760 μ L and 1.71 mL at a time is ample, as is described here.

3.1. Preparation of the Transcription Mix (T1)

1. To prepare 760 μ L T1 Mix, put into a microcentrifuge tube on ice the following components, in order, with gentle mixing after each addition (to avoid problems with precipitation): 487 μ L water, 40 μ L 1 M HEPES, pH 7.6, 3 μ L 2 M $MgCl_2$, 20 μ L 100 mM spermidine, 10 μ L 1 M DTT, 100 μ L 10X NTPs, 100 μ L 10X Cap. Mix well (*see Note 8*).
2. Prepare aliquots (100 μ L) into pre-chilled microcentrifuge tubes on ice, freeze in nitrogen, and store at -80°C . Do not freeze-thaw more than four times. Stable for at least 2 years.

3.2. Preparation of the Translation Mix (T2)

Before using the RRL, it must be supplemented with hemin [to prevent translational inhibition due to eIF2- α phosphorylation by the heme-regulated kinase (13)] and treated with micrococcal nuclease (to digest endogenous mRNAs, primary coding for globin). This is done on 10 mL at a time and sub-aliquoted for later use to make the T2. Steps 1–5 below describe how to nucleate the RRL, and steps 7–9 describe the preparation of T2.

1. Thaw a 10 mL aliquot of crude reticulocyte lysate quickly and put immediately on ice.
2. Add 400 μ L of hemin solution (final concentration will be 4 μ M), 100 μ L of 100 mM $CaCl_2$, and 100 μ L of micrococcal nuclease (15 U/ μ L stock). Mix gently but thoroughly (by repeated inversion).
3. Incubate in 25°C water bath for 12 min, making sure the entire sample is immersed in the water to ensure even warming. Mix gently by inversion after ~ 3 or 4 min of incubation.

4. Transfer to ice, immediately add 100 μ L of 200 mM EGTA, and mix gently but thoroughly by repeated inversion.
5. Dispense 1 mL aliquots into pre-chilled 2 mL *microcentrifuge tubes* on ice, thereby leaving enough room to add other components to prepare the T2 mix below.
6. Freeze in liquid nitrogen and store at -80°C . Nucleated RRL should generally be stable for at least 1–2 years.
7. To prepare T2 mix, thaw quickly and put immediately on ice a 1 mL aliquot of hemin/nuclease-treated RRL from step 6 above.
8. Add the following, in order: 6 μ L creatine kinase (20 mg/ mL), 30 μ L tRNA (10 mg/mL), 224 μ L water, 300 μ L Emix, and 150 μ L CB 20X. The total volume will be 1.71 mL. Mix gently but thoroughly (*see Note 8*).
9. Dispense 200 μ L aliquots into pre-chilled tubes on ice, freeze in nitrogen, and store at -80°C . The T2 mix is stable for up to a year at -80°C and will tolerate around four or five freeze-thaws if handled properly.

3.3. Linked Transcription and Translation Reactions

The translation mixture is optimized such that the products of the transcription reaction should be used directly in the translation reaction without any purification of the transcript. This ‘linked’ system works because the composition of the translation mix is adjusted to account for the Mg^{+2} , DTT, and spermidine contributed from the transcription. Thus, it is critical to use the transcription reaction directly (i.e., not purified transcript or mRNA) because it is providing Mg^{+2} , DTT, and spermidine to the translation, all of which are important. If purified RNA is used, then you need to add Mg^{+2} to 1.2 mM, spermidine to 0.4 mM, and DTT to at least 0.5 mM (to maintain a reducing environment in the translation reaction). The protocol below is an example of a typical linked transcription–translation reaction. It can be scaled as necessary.

1. Design and obtain oligos to PCR amplify the coding region of interest (typically from a plasmid-containing cDNA of interest). The PCR product should contain a 5' T7 promoter for transcription. The 3' primer should anneal at or beyond the stop codon. The following 5' oligo contains the T7 promoter (*italics*) followed by a few linker nucleotides, a Kozak's site, the start codon (**bold**), and the nucleotides that should be chosen to anneal to the coding region to be amplified (indicated here by underlined Ns) 5'-*TAATACGACTCACTATAGGGGAGACCATG*
NNNNNNNNNNNNNNNNNNNN-3'.
2. Use the above oligos and desired template plasmid to PCR amplify the coding region of interest. Any of several

thermostable polymerases can be used, following the manufacturer's supplied reagents and protocol. As an example, a typical 100 μ L PCR reaction for Taq polymerase (New England Biolabs) contains 1X PCR buffer (final concentration; supplied as a 10X stock with the polymerase), 20 ng template plasmid, 1 μ M (final concentration) of each primer, 200 μ M (final concentration) each dNTP, 5 units Taq polymerase, water to 100 μ L. We usually perform 30 cycles.

3. Use the Qiagen PCR purification kit to purify the PCR product as described by the manufacturer. Elute the product from the spin column with 50 μ L water (not TE). Check an aliquot (usually 1 or 2 μ L) on an agarose gel to confirm and estimate the quantity of amplification. Typically, 50 ng/ μ L concentration is normal for an average-sized product (\sim 750 bp). Anywhere from 25 to 200 ng/ μ L should suffice for use as template in the transcription reaction below (*see* **Notes 14** and **15**).
4. *For each DNA template to be tested, set up on ice the following:* 7.6 μ L T1 mix, 0.2 μ L RNAsin, 0.2 μ L RNA polymerase (either SP6 or T7, depending on the promoter of your template DNA), and 2 μ L template DNA (*see* **Note 16**).
5. Incubate for 60 min at 40°C for SP6 and at 37°C for T7, then transfer to ice (*see* **Note 17**).
6. *To translate the products of the transcription reaction, add the following directly to the completed transcription reaction on ice:* 28.5 μ L T2 mix, 5 μ L 35 S-Met, and water/other reagents as desired to a final reaction volume of 50 μ L. If co-translational translocation is being carried out, rough microsomes (RMs), proteoliposomes, liposomes, etc. are included, typically at \sim 2 to 5 μ L per 50 μ L reaction (*see* **Note 18**).
7. Incubate the translation reaction at 32°C for 30 min (for the typical 25–50 kD product). Longer for larger products, shorter for smaller products (*see* **Note 19**).
8. Following the translation reaction, transfer the tubes to ice. Remove 1 μ L to a separate tube containing 19 μ L SDS-PAGE sample buffer for direct analysis, keeping the remaining 49 μ L on ice for downstream analysis (*see* **Section 3.2.**). Typically, we run half of this (saving the other 10 μ L in case of technical problems with the gel) on a 0.75 mm thick minigel, which is then fixed, Coomassie stained (to confirm equal loading of all lanes), dried, and applied to film (Kodak-MR single emulsion film). Such direct analysis of the total translation products is helpful in troubleshooting downstream assays (*see* **Section 3.2.**) in case they yield confusing results because you will know exactly what you started with before any additional manipulations were performed (*see* **Note 20**).

3.4. Assays for Translocation and Topology

There are several ways to assay the segregation of proteins to the ER membrane. The three main ones are depicted in **Fig. 20.2**. Other assays have also been employed, but are not discussed further here. The first is to exploit an ER-specific modification such as glycosylation or signal sequence cleavage (**Fig. 20.2A**). Although this is very straightforward (typically detected as a

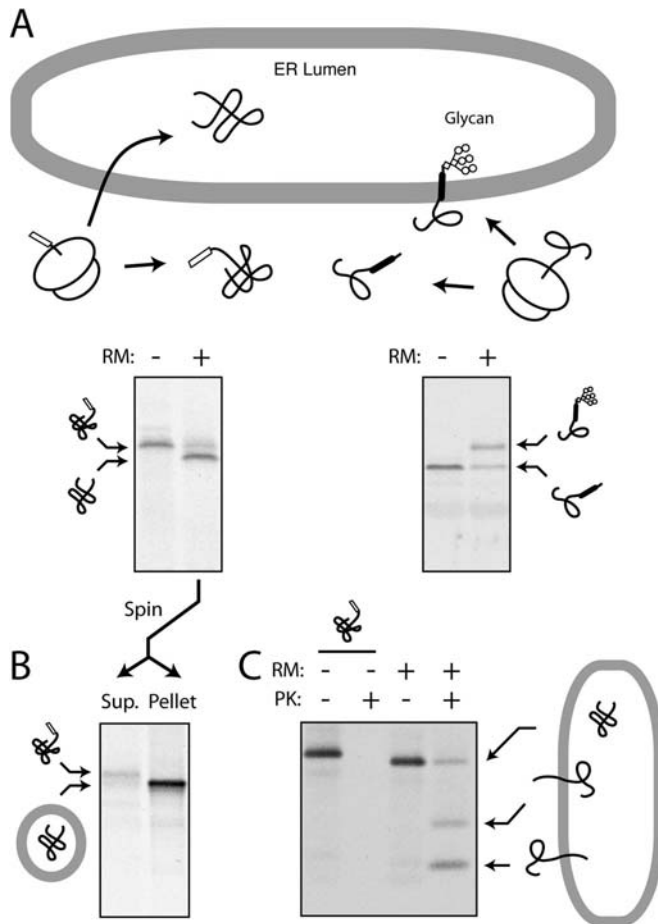


Fig. 20.2. Assays for protein segregation to the ER. **(A)** ER translocation-dependent modification. In the example on the *left*, a precursor becomes processed by signal peptidase only upon its translocation into the lumen of rough microsomes (RM), an event that can be monitored by a change in migration on SDS-PAGE. In the example on the *right*, a glycosylation site becomes modified upon successful translocation (8). **(B)** Co-fractionation assays. A sample similar to that from panel **A** can be separated by centrifugation into a cytosolic supernatant and membrane pellet to assess successful translocation. **(C)** Protease protection assay. Upon addition of proteinase K (PK) to the products of a translocation reaction, proteins that are either fully or partially translocated into the lumen of RMs are protected. Even a protein that generates multiple topological forms (such as mammalian prion protein; *see ref. 20*) can be resolved by this assay. By contrast, lack of translocation leads to complete digestion upon PK treatment.

change in migration by SDS-PAGE), it is not always applicable, may need the introduction of artificial glycosylation sites or epitopes into the substrate (8), and depends on the presence of functional and topologically restricted enzymes with very high activity and efficiency. This generally precludes its use if the ER components are fractionated. The second approach is to separate ER membranes from cytosol and determine whether the substrate cofractionates with the ER (**Fig. 20.2B**). Varying levels of stringency can be applied during the fractionation (e.g., high salt, pH 11.5, urea) to increase the specificity of substrate association and reduce background. However, even with these precautions, this approach cannot be used to definitively illustrate translocation (versus peripheral association) and cannot reliably provide information about protein topology. Nonetheless, such methods are very useful to re-isolate the membranes for post-translocational assays (e.g., retrotranslocation). And finally, protease protection assays can be employed to assay both translocation and topology (**Fig. 20.2C**). In this strategy, proteins (or portions of proteins) that are translocated or membrane inserted are shielded from proteases added to the cytosolic side of the membrane. By employing very high concentrations of an aggressive and relatively non-specific protease, essentially all cytosolically exposed protein can be digested to leave only the specific translocated population. By combining this approach with immunoprecipitation, the background can be markedly minimized to allow detection of even very low efficiency translocation events. Below is a typical protocol for both fractionation (by both sedimentation and floatation) and protease protection assays from the 50 μ L translation reaction produced in **Section 3.3**.

*3.4.1. Isolation of
Microsomes from
Translation Reactions by
Sedimentation*

1. Remove 10 μ L of the translation reaction and dilute with 90 μ L 1X PSB on ice.
2. Layer this onto 100 μ L of the sucrose cushion in an ultracentrifuge tube (either a tube for the TLA120.1 rotor or a micro-test tube for the TL100.3 rotor).
3. Spin for 5 min at 200,000 g (e.g., 70,000 rpms in the TL100.3 rotor). This spin time is suitable for traditional pancreatic RMs. If semi-permeabilized cells are used as source of ER, 5 min in a microcentrifuge may be sufficient. If proteoliposomes or smooth ER is used, longer spin times may be needed (e.g., 200,000 g for 30 min).
4. Remove the supernatant to a separate tube and resuspend the pellet in 1X PSB. It can be analyzed further if desired (e.g., by a protease protection assay as in **Section 3.4.3**.) or directly prepared for SDS-PAGE.

5. Analysis of the fractionation can be assessed by running equivalent amounts of the total (saved in step 8 of **Section 3.3.**), supernatant, and pellet fractions.

3.4.2. Isolation of Microsomes from Translation Reactions by Floatation

1. Remove 10 μL of the translation reaction and dilute with 90 μL of 2.2 M sucrose in 1X PSB. Considerable mixing will be required to ensure homogeneity.
2. Put into the bottom of a TLA120.1 tube. Layer with 100 μL of 1.8 M sucrose in 1X PSB followed by 25 μL of 1X PSB.
3. Centrifuge at 350,000 g for 1 h at 4°C.
4. Carefully remove the top 60 μL (the membranes will have floated to the top of the 1.8 M sucrose step). If the vesicles in this sample need to be recovered, the sample can be diluted in 1X PSB to 200 μL and centrifuged at 200,000 g for 20 min.
5. The floated vesicles can be analyzed by SDS-PAGE relative to the starting sample to assess the extent of membrane association.

3.4.3. Protease Protection Assay

1. Divide the translation reaction into three aliquots of 9 μL each on ice (A, B, and C).
2. Add 1 μL PSB to A, 0.5 μL PK and 0.5 μL PSB to B, and 0.5 μL PK and 0.5 μL Triton X-100 to C. Mix well and incubate on ice for 60 min.
3. To terminate the protease digest, both a protease inhibitor and rapid transfer to denaturing conditions are used. Either step alone is not fully sufficient to avoid artifacts. Approximately 10 min before the digestion reaction above is completed, start a boiling water bath.
4. Aliquot 100 μL of PKB into an appropriate number of empty tubes corresponding to the proteolysis reactions.
5. Dissolve a small amount of PMSF (2–5 mg) in DMSO to 0.25 M at room temperature.
6. Dip a P-2 pipette to draw up a very small amount (~ 0.1 – 0.2 μL) into the tip by capillary action and expel this into each proteolysis tube, mix, and put on ice.
7. After all tubes are completed, put the tubes containing the PKB into boiling water bath to pre-heat them to 100°C.
8. After tubes heated up (~ 1 – 2 min), transfer each proteolysis reaction directly into the boiling SDS solution and mix by rapidly pipetting up/down several times. Continue

boiling sample for additional 1–2 min and remove to room temperature.

9. Analyze an aliquot (5 or 10 μL) of the three reactions (A–C) by SDS-PAGE and autoradiography. The remainder can be subjected to immunoprecipitation if desired.

3.5. Reconstitution of Membrane Proteins into Proteoliposomes

The basic steps in preparing reconstituted proteoliposomes are to prepare a mixture of membrane proteins, detergent, and lipid. The detergent is then slowly removed, during which the detergent–protein–lipid micelles will assemble into lipid vesicles containing the membrane proteins (**Fig. 20.3**). Removal of the detergent can be accomplished in several ways including sim-

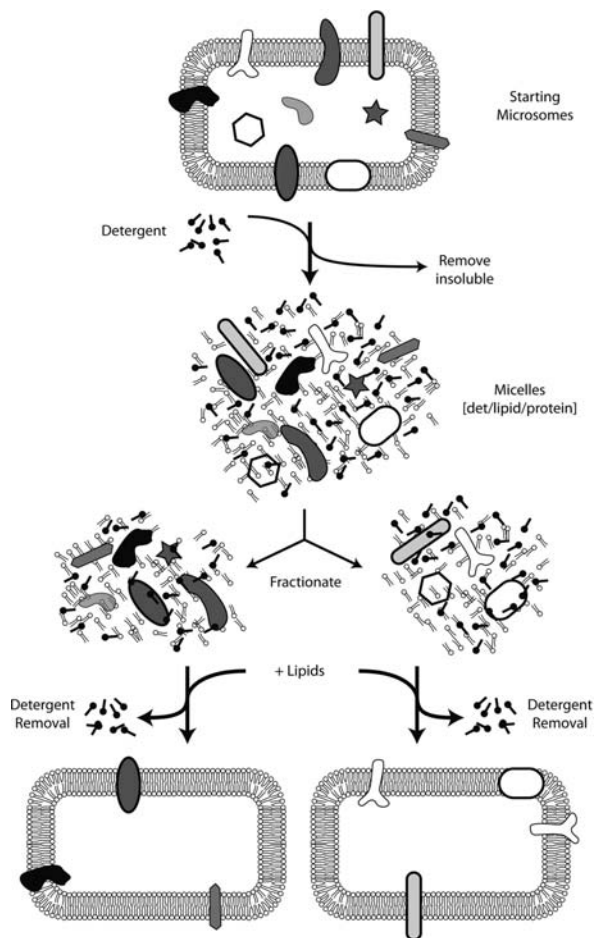


Fig. 20.3. Schematic depiction of membrane protein reconstitution. Crude microsomes are solubilized with detergent, fractionated, and reconstituted into proteoliposomes by removal of detergent in the presence of phospholipids. Note that not all proteins are successfully reconstituted, and the orientation achieved after reconstitution must be checked empirically.

ple dilution, dialysis, or adsorption. The method below is based on adsorption of the detergent to hydrophobic beads (Bio-Beads SM2, from Bio-Rad). The optimal conditions for reconstitution are difficult to predict and must be determined on a case-by-case basis. Furthermore, the orientation acquired by the membrane protein in the proteoliposomes is also not stochastic. Although it is often assumed that the orientation will be random, direct analysis shows this is not the case for many individual proteins.

From a practical standpoint, the protein, detergent, and lipids come from different sources and are mixed together just before reconstitution. The membrane protein(s), either a crude mixture or purified proteins, are already in detergent solutions to make them soluble. In addition, crude detergent extracts of membranes will also contain lipids. If the detergent and lipid contents need to be changed, the easiest method is to bind the membrane proteins to a chromatography resin, wash extensively with buffers containing the desired detergent, and elute in a buffer containing this new detergent. This is typically known as ‘detergent exchange’ and is often used if the detergent used for solubilization or purification is different from the one that proves to be best for reconstitution. A common example is the use of digitonin for solubilization and purification (due to its very gentle properties in maintaining membrane protein complexes), but exchange to another detergent just prior to reconstitution (because digitonin is very difficult to remove by dialysis or binding to Bio-Beads). *See ref. (11) for an example.*

The lipids are often provided separately, allowing the investigator to control the composition of the resulting proteoliposomes. They are usually easiest to prepare and handle as liposomes, but need to be added to the membrane protein as detergent-solubilized micelles. This way, the membrane proteins, lipid, and detergent will form mixed micelles at the start of the reconstitution.

3.5.1. Preparation of Lipids and Bio-Beads

Lipids are essential to making proteoliposomes containing solubilized membrane proteins. Although the precise composition of lipids may be largely irrelevant for reconstitution per se, they can significantly influence the activities of the reconstituted proteins. A simple mixture of phosphatidylcholine (PC) and phosphatidylethanolamine (PE) from a natural source (liver) in a 4:1 ratio is used here. However, other lipids, as well as cholesterol, can also be included if these are deemed important. Lipids are prone to oxidation and are therefore typically supplied in sealed ampules containing an inert gas (such as argon). With further manipulations, a reducing agent (e.g., DTT) is often included to prevent oxidation. In addition, lipids in organic solvents are typically not handled with plastic pipettes or put into plastic tubes. This is

because the lipids often stick to the plastic, and the organic solvent can extract various contaminants from the plastic into your sample. Once lipids are hydrated in aqueous solutions (in which they form liposomes) or solubilized in detergent (in which they form micelles), they can be handled with the usual plastic pipettes and tubes.

1. Using glass measuring tools (e.g., glass pipettes), transfer 20 mg PC and 5 mg PE to a glass test tube or glass vial. If a tracer is desired to follow lipid recovery, include 0.5 mg rhodamine-DPPE (with accordingly less PE).
2. Dry the lipids under high vacuum. A conventional SpeedVac can be used. We usually leave it overnight to ensure complete removal of solvent. Alternatively, the organic solvent can be evaporated and the lipids dried to a film on the side of the tube using a stream of nitrogen (do this in the hood to avoid breathing chloroform vapors). After the bulk solvent is removed, the lipids can be lyophilized under high vacuum to remove any traces of organic solvent.
3. Hydrate the lipids by adding 0.4 mL lipid hydration buffer to the dried lipids and resuspend by agitation, vortexing, and/or sonication in a bath sonicator. Resist the temptation to mechanically resuspend the lipids (e.g., with a pipette) as they can be quite sticky at this stage. Be patient, and they will be fully resuspended with vortexing, at which point you will have a homogeneous cloudy/milky suspension of liposomes. Once resuspended, they can be handled with plastic pipettes and put into plastic tubes.
4. Transfer the lipids to an Eppendorf tube and adjust to 0.5 mL with additional lipid hydration buffer as necessary to make a 50 mg/mL suspension. This can be divided into 100 μ L aliquots, frozen in liquid nitrogen, and stored at -80°C (*see Note 22*).
5. Put ~ 20 – 30 mL of dry Bio-Beads into a 50 mL polypropylene tube. Fill with MeOH, mix, and let the beads settle. Pour off the MeOH and fill with distilled water. Mix, let the beads settle (or if you are impatient, brief centrifugation), and pour off the water. Repeat this extensively (20 or more times) until all traces of MeOH are removed. Alternatively, pour the beads out into a disposable filter flask and use a vacuum to extensively wash the beads with distilled water. Put the beads back into a 50 mL polypropylene tube, fill with water, and store at 4°C . They are stable indefinitely.

3.5.2. Solubilization of ER Membranes

The solubilization of membrane proteins with detergent is influenced by many parameters. The most important are (i) the choice of detergent (*see Note 21*), (ii) the detergent concentration (more specifically, the relative ratio of detergent:protein:lipid),

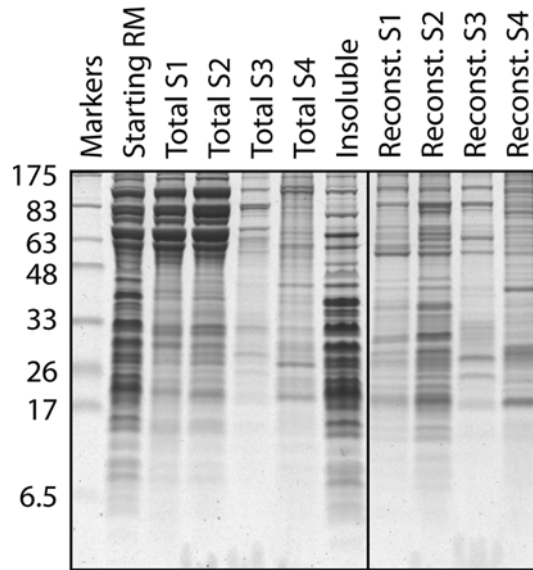


Fig. 20.4. Example of differential solubilization and reconstitution of membrane proteins. Crude RMs (*lane 1*) were sequentially extracted by four buffers containing different amounts of detergent and salt to generate four supernatant fractions (S1–S4) and an insoluble pellet (primarily containing ribosomes). Each of these four fractions was then reconstituted in the presence of phospholipids by detergent removal, and the resulting proteoliposomes were also analyzed on the gel. Note that the very abundant high molecular weight proteins in S1 and S2 (primarily luminal proteins) are not reconstituted. Note also the different protein profiles of the different proteoliposomes illustrating the utility of differential membrane protein extraction as a purification step. The detergent in this case was DeoxyBigCHAP, although similar results can be obtained with other detergents.

and (iii) the concentration and type of salt present during solubilization. The fact that different proteins are solubilized under different conditions can actually be used to effect some degree of purification via sequential selective solubilization (*see Fig. 20.4*). In addition to solubilization of proteins, maintaining its functionality further constrains the conditions. Hence, conditions that most efficiently solubilize a membrane protein may also result in its irreversible denaturation. While there are general guidelines, the solubilization conditions for any given protein or activity must be determined empirically. A general protein stabilizing agent such as 10–15% glycerol is often helpful. Some of the other less common parameters that influence stability of certain types of membrane proteins include divalent cations (e.g., Mg^{+2}), nucleotides, free phospholipids, or specific co-factors. Below is a generic protocol for solubilizing and enriching for most ER membrane proteins under conditions that allow their relatively straightforward reconstitution. Although not described, the solubilized proteins can of course be fractionated (e.g., by ion exchange or lectin chromatography) prior to reconstitution of individual fractions by similar methods.

1. 1 mL pancreatic rough microsomes (RM) at a concentration of 50 A₂₆₀ units is put into a 3 mL thickwall polycarbonate centrifuge tube (Beckman) on ice. Add an equal volume of ice-cold pre-extraction buffer and mix well.
2. Centrifuge at 400,000 *g* in a TL100.3 rotor for 15 min at 4°C. Remove the supernatant (which will contain primarily luminal proteins and some lipids that have been extracted by the low concentration of detergent making holes in the membrane). This can be saved for other uses if desired.
3. Resuspend the pellet in 0.9 mL ice-cold extraction buffer. This can be done by either gentle repeated pipetting or transfer to a small smooth glass homogenizer and manual homogenization with either a teflon or glass pestle. After resuspension, the solution should be homogeneous and turbid.
4. Add 100 µL of the 10% DBC on ice to solubilize the membrane. The solution should become more clear. Let sit on ice for 5–15 min to allow complete solubilization.
5. Centrifuge at 400,000 *g* in a TL100.3 rotor for 30 min at 4°C. Transfer the supernatant (the detergent extract, containing most of the membrane proteins and lipids) to a fresh tube on ice for use in the reconstitutions (*see* step 3 of **Section 3.3.**). The pellet will contain ribosomes, any tightly associated membrane proteins, and other membrane proteins that were not solubilized in step 4. If these remaining membrane proteins are desired, they can be solubilized in a subsequent step using harsher conditions (e.g., a different detergent, higher salt). Otherwise, it can be discarded.

3.5.3. Reconstitution of Solubilized Membrane Proteins

The amount of Bio-Beads needed to effectively remove the detergent without significant removal of the lipids or protein, while accomplishing this slowly enough to allow formation of proteoliposomes without aggregation, needs to be determined empirically. The optimum varies depending on the detergent being used, its concentration, and the salt conditions. Below is a sample protocol for how we usually do this for any given detergent extract or sample (*see* **Note 23**).

1. Mix 40 µL of the liposome stock (prepared in step 1 of **Section 3.3.**) with 10 µL of the 10% DBC. If your samples to be reconstituted were prepared in another detergent, use that to solubilize the liposomes. Thus, the sample to be reconstituted typically contains only a single detergent.
2. Add the solubilized lipids from step 1 to the 1 mL of detergent extract (prepared in step 2 of **Section 3.3.**) on ice.

3. Dispense different amounts of the Bio-Beads (prepared in step 1 of **Section 3.3.**) into 0.5 mL microcentrifuge tubes and remove the water. A good range to start is between 30 and 200 μL of packed bead volume. An alternative is to use round-bottomed 2 mL microcentrifuge tubes.
4. Add 105 μL of the detergent–lipid mixture (from step 2) to each of the tubes containing the Bio-Beads.
5. Incubate with gentle overhead mixing for ~ 12 – 16 h at 4°C . If using the 2 mL tubes, use orbital shaking to mix (e.g., an Eppendorf Thermomixer place in the cold room, set at 700–800 rmps).
6. Separate the fluid from the beads. If using the 0.5 mL tubes, we briefly spin down the sample (a pulse in the microcentrifuge), cut off the lid, puncture the bottom with a 26-gauge needle, put the tube into a larger 1.5 mL microcentrifuge tube, and spin for ~ 1 min to recover the fluid into the larger tube (while retaining the beads in the smaller tubes). Alternatively, remove the fluid carefully with a thin long pipette tip (e.g., a gel-loading tip) and transfer to 1.5 mL tubes.
7. Add five volumes (500 μL) of ice-cold water (or 1X PSB). This serves to dilute any residual detergent and reduce the concentration of glycerol to reduce the density of the solution.
8. Centrifuge the samples in a TL100.3 rotor with micro-test tubes and adaptors at 200,000 g for 30 min.
9. Remove the supernatant and resuspend the pellet in 25 μL 1X PSB by careful repeated pipetting.
10. The efficiency of lipid recovery can be monitored by including rhodamine-DPPE in the lipid mixture (*see* step 1 of **Section 3.3.**). The amount of lipid in the starting sample (step 2 above) can then be compared to the amount in the final proteoliposomes (step 9 above) by measuring absorbance at 560 nm.
11. The efficiency of protein reconstitution can be assessed by comparing the starting extract and final proteoliposomes by SDS-PAGE and staining for total proteins or immunoblotting for individual membrane proteins.

4. Notes

1. We use several precautions for handling all of the reagents for the T1 and T2 mixes. Use de-ionized, clean,

RNAse-free water for all of the solutions. In general (unless otherwise indicated), we quick-freeze items by immersion in liquid nitrogen and thaw quickly using either the warmth of your hands or a room temperature water bath. The rapid freezing and thawing are important for a few reasons. First, it minimizes/avoids precipitation of various reagents that can occur upon slow freezing/thawing. Second, it minimizes oxidation. Third, proteins are better protected from damage/denaturation by ice crystals if freezing/thawing is rapid. Immediately after thawing, reagents are generally kept on ice unless otherwise noted and immediately put away after use.

2. The pH of solutions is important. To avoid differences in the ways different people measure pH (different brands of pH paper, pH meter, etc.), we have provided the exact amounts of acid/base to add for most solutions that require accurate buffer/pH conditions. Other items can be adjusted to approximately the indicated pH using pH paper as a rough guide. We also typically favor HEPES over Tris buffers because the former has a pK_a closer to 7, and its pH does not vary with temperature.
3. Cap can be omitted from the transcription reaction if cost is an issue, and if an approximately twofold lower efficiency of translation can be tolerated. In this case, use 5 mM instead of 1 mM GTP in the 10X NTPs stock above.
4. Supplementing RRL translations with tRNA is not absolutely necessary and can be omitted if cost is an issue. It does, however, stimulate the translation of some proteins (presumably because the endogenous levels of some tRNAs in the RRL are limiting), and it is possible that some RRL preparations contain relatively low amounts of tRNA and thus would benefit from supplementation.
5. Endogenous RRL usually contains reasonable amounts of all of the amino acids, making supplementation optional. However, some batches of RRL may contain more or less of particular amino acids, so we always supplement with a complete mixture to avoid variability.
6. The most active and well-characterized source of ER-derived microsomes is from canine pancreas. Preparation of canine pancreas RMs in-house follows a well-established protocol (1). Microsomes have also been prepared from many other sources including rat liver (2) and cultured cells. In addition, semi-permeabilized cells have also been used and are a viable alternative (3). However, for the purposes of fractionation and reconstitution studies, rather large amounts of microsomes are needed and can realisti-

cally only be isolated from tissue since such large-scale cell culture is prohibitive. It is worth noting that microsomes from different sources may have different functional properties even though they are all ER-derived.

7. PK should be of the highest purity available commercially. Because it is purified from natural sources, contamination with other enzymes that could affect reliability or disrupt membranes (e.g., lipases) is a potential concern. Historically, preparations of PK were 'pre-digested' for 10 min at 25°C to proteolytically destroy any possible contaminants (PK is stable to its auto-digestion), but this has proven unnecessary in our experience.
8. Mix solutions gently but thoroughly. Keep in mind that many components are of considerably different density and do not always mix together with just a little tap of the tube. When mixing protein-containing solutions, avoid making bubbles or frothing the sample to minimize the possibility of denaturation.
9. RNAses are a potential problem because they are ubiquitous and will obviously preclude the transcription and translation reactions. We have found that simple cleanliness and care are the best strategies. The use of DEPC-treated water for the various reagents is helpful. Keeping your pipettes and tips clean and free from dust is helpful. Do not touch your tubes and tips with bare hands. The inclusion of RNasin is helpful, but only inhibits certain RNAses, and is not really necessary. In general, we can do all of this without RNase inhibitor by following the above precautions, and we do not have problems with degradation.
10. Optimal conditions for translation can vary somewhat between different transcripts. The most important parameters are Mg^{+2} (typically between 1.5 and 3 mM) and K^{+} (typically between 50 and 150 mM). Note that because a crude reticulocyte lysate is used, it contains endogenous Mg^{+2} and K^{+} (estimated to be ~ 1.7 and ~ 40 mM, respectively), probably at around one-third the levels typical for the cytosol. Hence the need to supplement these salts, especially because the added nucleotides can chelate Mg^{+2} .
11. Microsomes subjected to multiple freeze-thaws will lose luminal contents by leakage. Presumably, the integrity of the vesicles is sensitive to freeze-thawing. If the presence of these proteins in the reaction might pose a problem for your assays, it is helpful to re-sediment the microsomes and resuspend them in PSB just before use.

12. Cheaper alternatives to purified liver phospholipids are egg or soy-derived PC and PE. However, the acyl chain compositions can be quite different and may influence your proteins and/or activities. This needs to be tested empirically. In addition synthetic lipids of defined acyl chain composition are also available (e.g., DOPC, POPC) if you wish to precisely control lipid composition of the membrane. Other fluorescent lipids are also available to suit your requirements.
13. DBC is a detergent of the bile salt family. Others in the same family are Cholate, CHAPS, and BigChap. Among these, Cholate is the least expensive, but is also anionic and would therefore interfere with ion exchange fractionation. CHAPS is a good alternative that is relatively inexpensive and zwitterionic, making it more suitable for use in fractionation. Another commonly used detergent for reconstitution is octyl-glucoside because it is uncharged and is easily removed by either dialysis or Bio-Beads. A more thorough discussion of detergents is presented in ref. (14).
14. Although circular or linearized plasmids can be used for transcription/translation, PCR products are usually easiest and work well. In either case, you will need an SP6 or T7 promoter and a consensus Kozak's sequence at the start codon. A poly-A tail is not necessary for efficient *in vitro* translation. Even if there is no promoter in the plasmid containing the coding region of interest, you can simply encode the SP6 or T7 promoter and Kozak's sequence into the 5' primer used for PCR (as is described in step 1 of **Section 3.3.**). Note that the 5' primer can be designed to anneal internally to generate an N-terminally truncated protein. Similarly, 3' primers containing a stop codon can be designed within the coding region to generate a C-terminally truncated protein. A DNA template lacking a stop codon will result in a translation product that remains tethered to the ribosome via the last amino-acyl tRNA (15). Such truncated translation intermediates have many applications (for example, *see* ref. 16–19), but are not considered further here.
15. If you have numerous constructs in a vector that already has an SP6 or T7 promoter (as we do), you can use the same primers to amplify them all for translation: just use a 5' primer that anneals to the SP6 or T7 site and a 3' primer that anneals to the vector sequence beyond the open reading frame. The circular plasmids can also be used directly in the transcription reaction. For this purpose, use DNA that is free of RNase contamination and is at between 100

and 1000 ng/ μ L concentration (as is typical for standard Qiagen minipreps of high-copy plasmids).

16. In instances where many transcription reactions are being performed simultaneously, a master mix lacking the DNA can be prepared, aliquoted, and supplemented with the different DNAs. This can easily be scaled up or down as needed. A control reaction with water instead of DNA can be performed to determine the extent of background in the translation reactions.
17. Transcription reactions can be incubated for longer times (up to 2 or 3 h), but there is no significant increase in yield under these conditions (the nucleotides become limiting and the pyrophosphate that is generated inhibits the polymerase). After transcription, the sample can be kept on ice for some time (an hour or two) without any concern of degradation. If frozen in nitrogen and stored at -80°C , the transcript can be used again. However, we almost always make it fresh for each experiment and discard any leftover material (to avoid the risk of degradation with storage).
18. The amount of ^{35}S -Met can be decreased depending on the amount of labeling needed. If unlabeled translations are being performed, cold Met (at a stock concentration of 0.4 mM) should be added instead (to 40 μM final concentration).
19. Incubations are generally at $\sim 30\text{--}32^{\circ}\text{C}$, but translation works at anywhere from 23 to 37°C . Typical translation times are 30–60 min, depending on size of expected products. In general, it takes ~ 5 min to complete a 100 residue protein (~ 10 kD) at 32°C . Thus, we generally incubate for ~ 10 min per 10 kD of expected size product. Longer incubations generally do not produce more product since translational activity of the lysate declines. Translation is slower at lower temperatures.
20. We find that for a small- to average-sized protein ($\sim 25\text{--}30$ kD) containing an average proportion of methionines, a band corresponding to the translation product can be detected in as little as 5 or 10 min for optimally translating proteins. Preprolactin and prion protein are often used in our lab, and both express comparably well. Note that hemoglobin, present in the RRL at ~ 50 mg/mL, migrates at 14 kD. Although it is not radio-labeled, it distorts this region of the gel and can cause artifacts if your translated product is of the same size. This can be avoided by loading 10-fold less of the translation product and exposing the gel for longer times. Alternatively, the translation products can be immunoprecipitated.

21. There have been many instances where detergents vary sufficiently from lot to lot to affect the results of functional reconstitution (e.g., ref. 11). It is worth keeping this in mind if you have issues with reproducibility. Some detergents are more variable than others, and ones from natural sources (e.g., digitonin) are especially difficult in terms of reliability. Cost can also be an issue in detergent choice. Thus, some degree of screening can be important in initial studies. The time invested is usually very worthwhile and pays off in the long run.
22. Although freezing and thawing liposomes are not recommended in many protocols, this is not really relevant here (as long as oxidation is avoided), since they will be subsequently solubilized in detergent. Thus, alteration of liposome size or maintaining them as unilamellar vesicles is not an issue.
23. Reconstitution of purified membrane proteins in detergent solution can be accomplished using basically the same methodology. However, it may be necessary to re-optimize conditions to maximize recovery, which can be expected to be at least 50%. Generally speaking, one tries to use the least amount of lipids needed to fully incorporate the protein of interest into vesicles (to maximize density in the membrane).

Acknowledgments

Work in the Hegde lab is supported by the NICHD Intramural Research Program of the National Institutes of Health.

References

1. Walter, P., and Blobel, G. (1983) Preparation of microsomal membranes for cotranslational protein translocation. *Methods Enzymol* **96**, 84–93.
2. Adelman, M.R., Blobel, G., and Sabatini, D.D. (1973) An improved cell fractionation procedure for the preparation of rat liver membrane-bound ribosomes. *J Cell Biol* **56**, 191–205.
3. Wilson, R., Allen, A.J., Oliver, J., Brookman, J.L., High, S., and Bulleid, N.J. (1995) The translocation, folding, assembly and redox-dependent degradation of secretory and membrane proteins in semi-permeabilized mammalian cells. *Biochem J* **307**, 679–687.
4. Stefanovic, S., and Hegde, R.S. (2007) Identification of a targeting factor for posttranslational membrane protein insertion into the ER. *Cell* **128**, 1147–1159.
5. Schuldiner, M., Metz, J., Schmid, V., Denic, V., Rakwalska, M., Schmitt, H.D., Schwappach, B., and Weissman, J.S. (2008) The GET complex mediates insertion of tail-anchored proteins into the ER membrane. *Cell* **134**, 634–645.
6. Daniels, R., Kurowski, B., Johnson, A.E., and Hebert, D.N. (2003) N-linked glycans direct the cotranslational folding pathway of influenza hemagglutinin. *Mol Cell* **11**, 79–90.

7. Oberdorf, J., and Skach, W.R. (2002) In vitro reconstitution of CFTR biogenesis and degradation. *Methods Mol Med* **70**, 295–310.
8. Brambillasca, S., Yabal, M., Soffientini, P., Stefanovic, S., Makarow, M., Hegde, R.S., and Borgese, N. (2005) Transmembrane topogenesis of a tail-anchored protein is modulated by membrane lipid composition. *EMBO J* **24**, 2533–2542.
9. Fons, R.D., Bogert, B.A., and Hegde, R.S. (2003) Substrate-specific function of the translocon-associated protein complex during translocation across the ER membrane. *J Cell Biol* **160**, 529–539.
10. Garrison, J.L., Kunkel, E.J., Hegde, R.S., and Taunton J. (2005) A substrate-specific inhibitor of protein translocation into the endoplasmic reticulum. *Nature* **436**, 285–289.
11. Görlich, D., and Rapoport, T.A. (1993) Protein translocation into proteoliposomes reconstituted from purified components of the endoplasmic reticulum membrane. *Cell* **75**, 615–630.
12. Jackson, R.J., and Hunt, T. (1983) Preparation and use of nuclease-treated rabbit reticulocyte lysates for the translation of eukaryotic messenger RNA. *Methods Enzymol* **96**, 50–74.
13. Trachsel, H., Ranu, R.S., and London, I.M. (1978) Regulation of protein synthesis in rabbit reticulocyte lysates: purification and characterization of heme-reversible translational inhibitor. *Proc Natl Acad Sci USA* **75**, 3654–3658.
14. Helenius, A., and Simons, K. (1975) Solubilization of membranes by detergents. *Biochim Biophys Acta* **415**, 29–79.
15. Perara, E., Rothman, R.E., and Lingappa, V.R. (1986) Uncoupling translocation from translation: implications for transport of proteins across membranes. *Science* **232**, 348–352.
16. Wiedmann, M., Kurzchalia, T.V., Hartmann, E., and Rapoport, T.A. (1987) A signal sequence receptor in the endoplasmic reticulum membrane. *Nature* **328**, 830–833.
17. Görlich, D., Hartmann, E., Prehn, S., and Rapoport, T.A. (1992) A protein of the endoplasmic reticulum involved early in polypeptide translocation. *Nature* **357**, 47–52.
18. Crowley, K.S., Reinhart, G.D., and Johnson, A.E. (1993) The signal sequence moves through a ribosomal tunnel into a noncytoplasmic aqueous environment at the ER membrane early in translocation. *Cell* **73**, 1101–1115.
19. Hegde, R.S., and Lingappa, V.R. (1996) Sequence-specific alteration of the ribosome-membrane junction exposes nascent secretory proteins to the cytosol. *Cell* **85**, 217–228.
20. Hegde, R.S., Mastrianni, J.A., Scott, M.R., DeFea, K.A., Tremblay, P., Torchia, M., DeArmond, S.J., Prusiner, S.B., and Lingappa, V.R. (1998) A transmembrane form of the prion protein in neurodegenerative disease. *Science* **279**, 827–834.